

Reference washer extractor Wascator FOM71 CLS

Features and benefits

- · Meets IEC and ISO international standards
- Suitable for testing washing effects of detergents and chemicals and for textile quality control
- Suspended drum construction allows a high extraction force.
 No foundation required
- Small space between outer/inner drum for efficient use of water and detergent
- · Dual water control of volume (weight) and level
- · Built in scale for precise volume control, accuracy:
 - volume control +/- 0.2 litre/bath
 - level control +/- 0.8 litre/bath
- Equipped with Clarus Control* a fully programmable electronic timer
- Frequency controlled motor system for flexibility in programming and precise speed
- · Test tap for water samples
- · Stainless steel drum assembly, front, side and top panel

Images shown are a representation of the product only and variations may occur.

Main options

- · WPM (Wash Program Manager) PC programming kit in Windows
- · Extra program memory cards

Main specifications			Wascator FOM71 CLS
Drum,	volume diameter	litre ø mm	61 520
		rpm	1100 350
Heating	electricity	kW	5.4



Function/subsystem


Ambient conditions	
The machine shall be able to function during the same ambient conditions as normal washing machines. The functional values with tolerances will be valid if following ambient conditions are kept	
Voltage Frequency Inlet water Cold Hot Water pressure	Nominal ±2% Nominal ±2% (15±2) °C Max (60±2) °C (240±50) kPa
Ambient temperature	(20±5) °C
Inner drum	4=
Diameter Depth Volume Lifting vanes: Number Height Radius Base width Perforation: Diameter	(520±1) mm (315±1) mm 61 l 3 53 mm 17 mm 65 mm 5 mm
Depth of countersink Total hole area Material: Stainless steel	2.5 mm 600 cm² 2333 (SS 18/8)
Outer drum	
Diameter Sump Material: Stainless steel	554 mm ~ (1300±25) ml 2333 (SS 18/8)
Door opening	
Opening diameter	310 mm
Cabinet	
Stainless steel	2333 (SS 18/8)
Drum speed	
Washing speed: Programmable Extraction: Programmable	(20 - 59±1) RPM (200 - 1100±20) RPM
Reversing rythm	
Programmable	(0 - 250±0.1) = sec. (Tolerance area refers to timer intervals)
Water level	
Weight sensing Repeatability Stepsize Dosing accuracy weight control	±0.1kg/bath 0.05kg ±0.2kg/bath
Level sensing: Repeatability Stepsize Dosing accuracy level control	±5 mm 3 mm ±0.8litre/bath
Thermostat	
Continuos variable in each wash sequence Selectable interval Programming steps Switch off temperature accuracy in Interval 30-97 °C Interval switch off temperature – switch on temperature	Yes (4-97) °C 1 °C ±1 °C 4 °C







Electrical o	connections	Wascator FOM71 CLS	
Heating	Voltage		
El.	230V 3AC 50 Hz	kW(A)	5.7(16)
	400-415V 3N AC 50		5.7(10)
415V 3N AC 50 Hz kW(5.7(10)
Water and	steam connection		
Water valves Water pressure Capacity at 300 kPa Drain valve Draining capacity Liquid detergent supplies		DN kPa l/min ø mm l/min	20 150-400 20 50/75 150 4
Floor requi	irements		
Frequency of the dynamic force Floor load at max extraction		Hz kN	18.3 1.6 ± 0.5
Sound leve	els		
Airborne sound level		dB(A)	71
Shipping d	ata		
Weight			195 260
Shipping vo	olume	crated, kg m ³	0.89
Accessorie	s		
Hose connection kit for water			X
Dimensions	s in mm		
A Width B Depth C Height D E F G H I K L M N O P R			720 690 1315 200 750 555 720 1230 220 1010 135 825 45 300 240
2 Cold wate 3 Hot wate 4 Drain		-	6 Control panel 7 Detergent compartment 8 Door opening Ø 310 mm


Other voltages available, see installation manual.

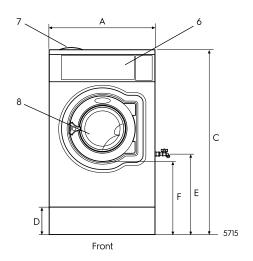

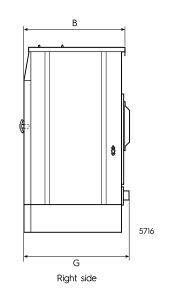

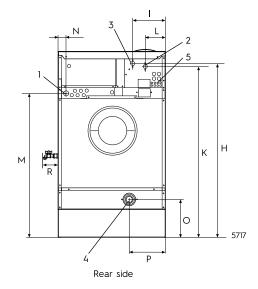



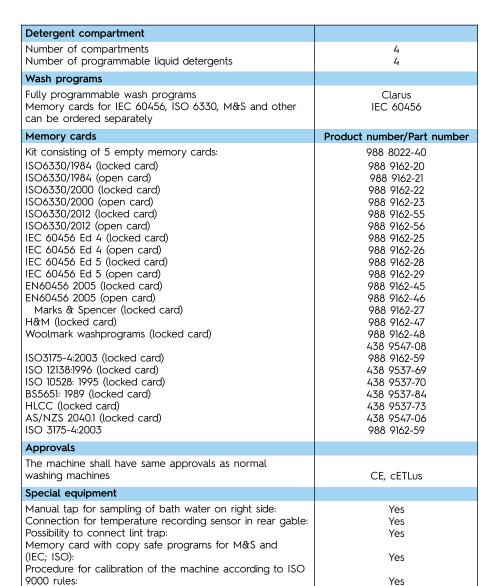



Approval from Marks & Spencer required.

wash:

systems:

Possibility to take out test pieces from wash bath during the

Analog output from load cells to data acquisition

Pre- programmed wash programs

The Clarus Control in the Wascator FOM71 CLS contains 9 pre- programmed "wash programs". They are not real wash programs but are mainly designated for checking purposes. The stored programs are:

purposes. The stored programs are:					
Place	Name	Content			
991	Start-up Programme	Flushes the machine with water. Shall be used if the machine has been idle for more than 2 hours.			
992	Level check	Checking of the level. With ruler inside the machine fills to 100 mm, 130 mm, 160 mm, and 200 mm.			
993	Weight check	Machine fills water to 18kg and to 26kg.			
994	Temperature check	Check of the temperatures at 25°C, 40°C, 60°C and 90°C 2 minutes at each temperature.			
995	Extraction	500 RPM 1 min.			
996	Extraction	500 RPM 5 min.			
997	Extraction	775 RPM 2 min.			
998	Extraction	950 RPM 4 min.			
999	Extraction	1100 RPM 5 min.			





Yes

Yes

Art No. 438918101/2015.01.14

the right to alter specifications without notice

Unit 3 / 41 Biscayne Way, Jandakot, WA 6164 Ph: (08) 9414 9474 Fax: (08) 6210 5961 Service: 0411 883 346 / 0400 246 352 Sales: 0406 659 429/ sales@comls.com.au

Web site : www.comls.com.au